


# EcoPac®

## Tip Refill System for Biomek® FX and NX Workstations

The benefits of pipette tip refills are now available to users of Biomek® FX and NX Workstations. Save on storage space, waste and costs with the new EcoPac® refill system from Alpha Laboratories.


The EcoPac system utilises environmentally friendly refillable racks and refill stacks of 5 x 96 robotic tips. They are compatible with instruments using the AP96 liquid handling head including

- Biomek® FX
- Biomek® NX
- Multimek®
- Biomek® 3000

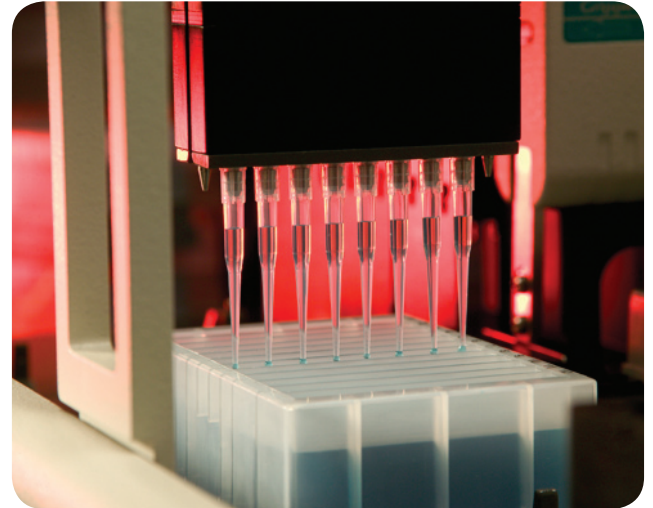

The refillable racks are supplied with a full tray of 96 tips and are used in place of the standard racked tips in your Biomek workstation. When the racks are empty, it is quick and easy to reload them using the EcoPac refill stack.

**Save Time** - It takes seconds to refill your racks

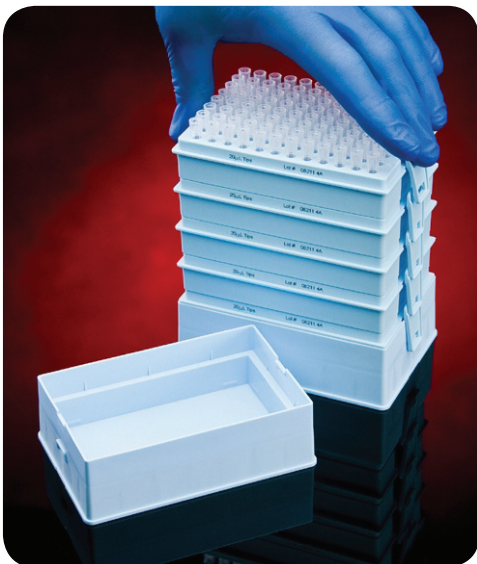

**Step 1:** Press in the tabs at each end of the empty tip tray and remove it from the rack.



**Step 2:** Lift a new tray of tips off the EcoPac refill stack using the same tab system to release it.



**Step 3:** Place the new tray into the empty rack and the tabs will lock it in position ready for use.

# EcoPac®

## Eco benefits and ordering details

### Save on Space and Waste

The compact EcoPac refills take up reduced storage space and generate less waste than traditional racked tips so are ideal for high usage laboratories.

### Save Money

Refills offer attractive cost savings compared to racked tips, allowing your budget to stretch further without compromising on quality.

### Peace of Mind

EcoPac tips are made from 100% virgin polypropylene and are metal-free. They are certified free from DNA, DNase, RNase and PCR inhibitors. Sterile tips are also tested for the presence of endotoxins.



### Compatibility Assured

EcoPac tips are subjected to rigorous physical testing. Each batch is checked for correct functionality using the original manufacturer's robots so you are assured of a trouble-free transition. There is no need to set up new tips in the instrument software. Just select your usual tip volume from the predefined list. See the ordering table below for further details.



## Ordering Details

Product Code	Description	Pack Size	Software Designation	Rack Colour Coding	Recommended maximum volume
LW1513R	20µl EcoPac Refillable Racked Tips	10 x 96	P20	Blue	70µl
LW1513T	20µl EcoPac Refills in 2 Stacks	10 x 96	P20	Blue	70µl
LW1513RS	20µl EcoPac Refillable Racked Tips, Sterile	10 x 96	P20	Blue	70µl
LW1513TS	20µl EcoPac Refills in 2 Stacks, Sterile	10 x 96	P20	Blue	70µl
LW1515R	250µl EcoPac Refillable Racked Tips	10 x 96	P250	Green	250µl
LW1515T	250µl EcoPac Refills in 2 Stacks	10 x 96	P250	Green	250µl
LW1515RS	250µl EcoPac Refillable Racked Tips, Sterile	10 x 96	P250	Green	250µl
LW1515TS	250µl EcoPac Refills in 2 Stacks, Sterile	10 x 96	P250	Green	250µl

### More Information

To discuss your requirements, test a sample or request a quotation, please contact us on 0800 387732 or e-mail [sales@alphalabs.co.uk](mailto:sales@alphalabs.co.uk).